


KaVo MASTERspace®

Prosthetic perfection,  
realized with  
KaVo MASTERspace!

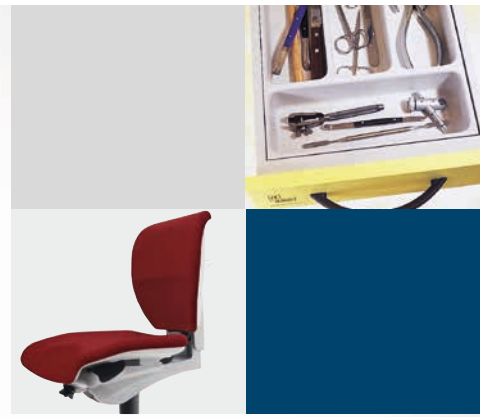

KaVo. Dental Excellence.

# Benefit from greater efficiency.

## Discover greater efficiency and motivation for your dental laboratory.

As one of the world's largest system providers, KaVo has accumulated decades of competence in high-quality laboratory equipment. With a high-performing team of experienced planning specialists, we develop installation solutions that are exactly customised to meet your requirements.

In this context increased efficiency and the motivation of your co-workers are central for our approach. We meet this high demand by optimising procedures as well as by ergonomic and functional planning and design of your workplaces.

As a laboratory technician, high-standards and precision come as second nature. Therefore the KaVo MASTERspace® is the natural choice.

You get your perfect laboratory equipment in the three equipment lines presented on the following pages:

KaVo MASTERspace® Select  
KaVo MASTERspace® Classic  
KaVo MASTERspace® Basic



# The three equipment lines.

KaVo offers you three equipment lines for your new laboratory installation. Depending on your requirements, decide in favour of MASTERspace Select, MASTERspace Classic or MASTERspace Basic.



- KaVo MASTERspace® Basic**
- Convincingly economical**
- Single-workplace solution without dust extraction.
  - Material cabinet including three powder-coated metal drawers.
  - Front selectable from standard colours in single-colour, handles matt chromium plated.



- KaVo MASTERspace® Classic**
- Souverän in Form und Funktion**
- With SMARTair® single-place dust extraction.
  - With EXTRAmatic® dust extraction unit for three workplaces.
  - Suction hood front table.
  - Front colours selectable from standard colours in single- or two-colour, handles matt chromium plated



- KaVo MASTERspace® Select**
- Perfect in function and design**
- With EXTRAmatic® multi-user dust extraction unit for the single and the three-person workplace.
  - Retractable armrests, suction hood and functional drawer
  - Front selectable.

Equipment lines	MASTERspace Basic	MASTERspace Classic	Select
SMARTair for single-place dust extraction	–	•	–
EXTRAmatic for single-place dust extraction	–	–	•
EXTRAmatic for multiple-place dust extraction	–	•	•
Retractable armrests	–	–	•
Retractable extraction hood	–	–	•
Extraction hood front table	–	•	–
Functional drawer	–	–	•
Front colours	Single-coloured	1-/2-coloured	Multicoloured
Handles	Metal handles chromium matt	Metal handles chromium matt	Metal handles chromium matt

• Standard version – not availableh

# KaVo MASTERspace Select.

## Timelessly beautiful design.

KaVo MASTERspace Select is the equipment line offering you perfection in function and design.

The functional drawer as well as the retractable armrests and the extraction hood are examples of highest functionality down to the smallest detail.

Both the single workplace and the three-person workplace of the MASTERspace Select are equipped with the multi-place dust extraction unit EXTRAmatic.

**MASTERspace Select:**  
Experience perfection.



Retractable armrests



Retractable extraction hood



Three-person workplace MASTERspace® Select withEXTRAmatic® and design package.



Single workplace MASTERspace® Select withEXTRAmatic® and design package.


# KaVo MASTERspace Classic.

## Convincing to the smallest detail.

KaVo MASTERspace® Classic is the economical alternative for all workplace solutions where dust extraction equipment is mandatory.

The single workplace is equipped with the SMARTair® single dust extraction unit, the three-person workplace with the EXTRAmatic® multiple-place dust extraction unit. MASTERspace® Classic is also available with function package and design package.

**MASTERspace Classic:**  
Work confidently.

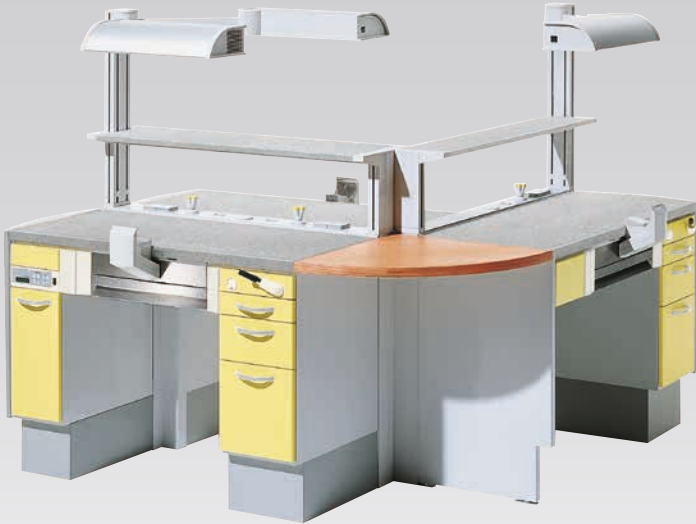

Three-person workplace MASTERspace® Classic withEXTRAmatic® and design package.



Single workplace MASTERspace® Classic with single dust extraction SMARTair® and function package.



# KaVo MASTERspace Basic.

## Working without dust extraction.

KaVo MASTERspace Basic is the ideal solution for single workplaces where dust extraction is not required. The equipment line Basic fits perfectly into the MASTERspace generation and can thus be ideally combined with KaVo MASTERspace Select or Classic. A material cabinet including three powder-coated metal drawers offers sufficient storage space. The waste drawers are made of high-quality stainless steel.

As a matter of course the Basic line is also available with function and design packages

## MASTERspace Basic: Convincingly economical.



MASTERspace® Basic single workplace with function package.



MASTERspace® Basic single workplace standard design.

# Convincing function and design.

With the KaVo function and design packages, you supplement your MASTERspace laboratory equipment quite specifically depending on workplace requirements and your personal style.

These optional packages can be used for all three lines, MASTERspace Select, Classic, and Basic.

An example of the enhanced functionality of the singular design is given by the KaVo workplace lamp

Due to its illuminance and contrast, the KaVo operating light is indispensable for all high-precision jobs.

- 2 x 36 W, Electronic ballast, 2,300 LUX performance (measured at a lamp-table distance of 62 cm).
- Distinctly longer service life by temperature-controlled fan cooling.
- Even light distribution and low reflection due to light grid.

## MASTERspace design packages



Design package single workplace consisting of single workplace lamp, worktray board, storage board, back screen and arch screens.



Design package group workplace consisting of single workplace lamp, double workplace lamp, 2 worktray boards, 3 storage boards and 3 arch screens.

## Highlights and advantages for you:

- Thanks to the individual KaVo 3-D laboratory planning, the work sequences are optimised, and the efficiency of the dental laboratory is increased.
- Superior design with the highest quality of workmanship and conservation of value.
- Optimum workflow is ensured by thought-out arrangement of the functional elements.
- Technically most advanced suction technology with KaVo SMARTair and KaVo EXTRAmatic for optimum functionality and clean laboratory air.

## MASTERspace function packages



Function package single workplace: Consisting of single-place lamp and work tray board.


Function package group workplace: Consisting of single-place lamp, two-place lamp and 2 work tray boards.

# Breathe safely.

## The lung is endangered.



The lung is particularly endangered by very fine dust.  
Feinstaub ist der Anteil des gesamten einatembaren Staubes im Labor, dessen Teilchen kleiner als 5 µm sind. Durch die Feinheit kann er mühelos bis in tiefste Luftwege, den Alveolarbereich, gelangen und dort als Folge ständiger Gewebereaktion die gefürchtete Staublunge auslösen.

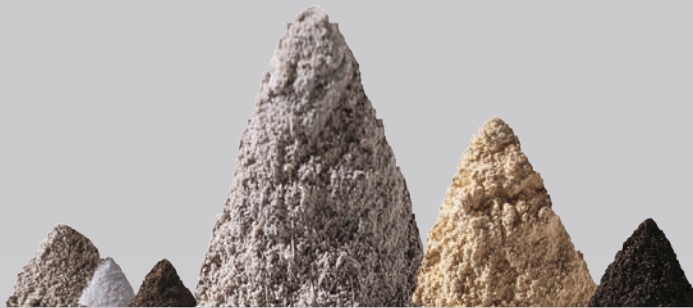

Ceramics steel noble metal methylmetacrylate gypsum others

Protect yourself and your laboratory personnel with KaVo dust extraction units. *Ø dust exposure/year and dental technician.*

KaVo laboratory extraction units are the result of over 20 years of experience in laboratory dust extraction – for your health.

## KaVo EXTRAmatic: The multiple-place dust extraction unit

Dust extraction is something many can do. However, the real benefit of a perfect dust extraction unit becomes evident only in the combination of optimum suction power and high functionality. Whether at the workplace proper or at the devices in the functional area does not matter.


- The ABSORBA long-term filtering system with integrated duplex safety filter.
- With the optionally inserted activated carbon filter, your KaVo EXTRAmatic® also simultaneously removes quite varied materials as well as plastic vapours.
- Thanks to the switch-on automatic, the suction motor works only when air is required and thus saves valuable energy as opposed to central extraction systems.
- The intelligent KaVo self-cleaning system saves filter and disposal costs day after day.
- If desired, purified exhaust air is fed back, thus saving money for heating your laboratory.
- The extraction level can be individually programmed for each of the up to three extraction locations and thereby offers perfect extraction management.
- Useful accessories enable integration of the KaVo EXTRAmatic® into any laboratory equipment program, any other workplace and best of all into the KaVo MASTERSpace® and FLEXspace® concepts.



## NEW KaVo SMARTair Evo: The single-workstation dust extraction unit



- Suction performance tested by DGUV
- Improved suction pressure and volume flow
- Dust separation with three-stage filter system
- Integrated automatic start and stop
- Automatic suction power adjustment
- Integrated motor protection
- Simple motor change
- Integration into third-party tables.

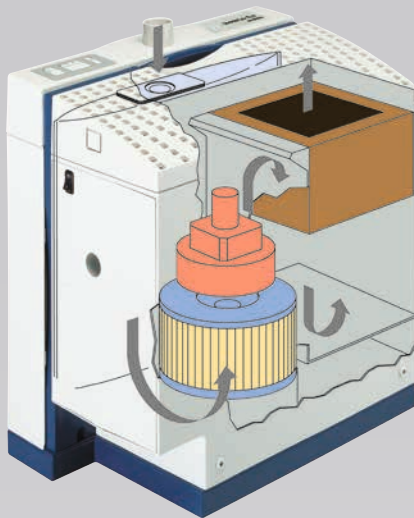

The integrated solution is ideal for technicians' tables.  
The drawer concept offers maximum comfort for the technician.



The mobile solution is primarily used in the functional area of the dental laboratory and also where spatial restrictions prevent an integrated solution.



# Sitting posture must be healthy.

## KaVo SENSit® – tailor-made for the dental laboratory.

The KaVo engineers have developed a chair that takes into account all health aspects of the sitting position. The anatomical design of the seat and the backrest supports the body precisely where it is necessary. And the shape and functionality of SENSit are ideally adapted to the working conditions in the laboratory. Ingeniously thought out to the smallest detail. For a relaxed, body-friendly sitting position while maintaining maximum freedom of movement.

**Material**

Slip-resistant, moisture-permeable, easy to clean, kind to the skin. The seat and backrest cushions can be very easily exchanged.

**Seat**

The seat precisely adapts to the individual anatomy of the person. The upholstery is specially designed for the typical bent posture working position of the dental technician. This avoids excessive pressure. You can move the seat forward and backward.

**Variable**

With the SENSit you always sit at the right height.

**Backrest shape**

The backrest shape is specially designed for work in the laboratory, it supports the spine and supports your back comfortably and safely.

**Height adjustment**

The backrest height can be adjusted so that the upholstery is always in the right area of the back.

**Synchronous mechanism**

Backrest and seat are articulately jointed. When the backrest opening angle increases, then the seat also bends back.

## KaVo SENSit® – the fine adjustment of the perfect seating position.

Have a seat! And do not hesitate to lean back and relax! Studies have shown: When sitting for a longer period of time, then changing the sitting posture is particularly important.

The seat can be moved horizontally so that the contact with the backrest is optimally adapted to the physique.

Change of position made easy: You determine to what extent seat and backrest bend! Lift handle. Strain backrest. Lock mechanical system. Done!

The height adjustment of the backrest permits individual adjustment of the chair to your body height.

Whether short or tall, the spine is always properly supported: Thanks to the height-adjustable backrest the special upholstery is where it should be.

You can individually adjust the resistance of the backrest: Turn the adjusting wheel and set the force of the synchronous mechanism to your body weight. Thus you always sit comfortably.

The robust aluminium pedestal: High load-bearing capacity and long service life due to excellent material quality.


# Here you have the choice.

## MASTERspace


### Table top coatings




ROCKdent



FORMdent  
Agate grey




Laminat  
Carrara

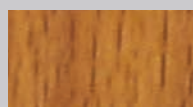

Multiplex  
Beech



Corian®  
dark

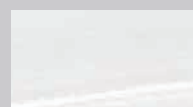


Corian®  
bright


### Fronts



Gentian blue



Tourmaline-anthracite



Carrara



Sunflower



Orange



Bordeaux

### Base




Stone grey

## SENSIT

### Fabric "Extreme"



stahlblau/  
steel blue

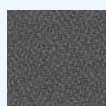

titangrau/  
titanium grey




korallrot/  
coral

- 100% Xtreme FR
- 400g/m


### Fabric "Microtouch"



saphirblau/  
sapphire

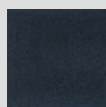

tiefseeblau/  
ocean blue

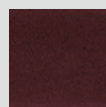

aubergine/  
aubergine

- 88% Polyester
- 12% Polyurethane/PU
- 385 g/m


### Fabric "Eccos"



olivingrün/  
olive green




polarsilber/  
polar silver

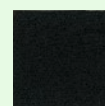

schwarz/  
black

- 9% Ramie
- 91% Wool
- 546 g/m

### "Leather"

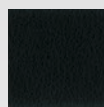

schwarz/  
black

- 100% Leder/  
leather



KaVo. Dental Excellence.